



Integra®
Cranial Stabilization
Product Catalogue

Table of Contents

Ordering and Warranty Information	4
MAYFIELD® Triad Skull Clamp	5
MAYFIELD® Base Units	6
MAYFIELD® complete Skull Clamp system	6
MAYFIELD® Ultra 360 Patient Positioning System	7
MAYFIELD® Ultra Base Unit	7
MAYFIELD® Swivel Adaptors	8
MAYFIELD® Swivel Adaptor	8
MAYFIELD® Modified Swivel Adaptor	8
MAYFIELD® Ball Socket Swivel Adaptor	8
MAYFIELD® Tri-Star Swivel Adaptor	9
MAYFIELD® Modified Tri-Star Swivel Adaptor	9
MAYFIELD® Skull Clamps	10
MAYFIELD® Infinity Skull Clamp system	10
MAYFIELD® Triad Skull Clamp	10
MAYFIELD® Modified Skull Clamp	11
MAYFIELD® 2000 Skull Clamp	11
MAYFIELD® Skull Pins	12
MAYFIELD® Reusable Skull Pins (Steel)	12
MAYFIELD® Disposable Skull Pins (Steel)	12
MAYFIELD® Disposable Skull Pins	12
MAYFIELD® Radiolucent Skull Pins	13
MAYFIELD® Disposable Titanium Skull Pins	13
MAYFIELD® Reusable Titanium Skull Pins	13
MAYFIELD® Horseshoe Headrests	14
MAYFIELD® complete Headrest System	14
MAYFIELD® Horseshoe Headrests	14, 15
MAYFIELD® Papadopoulos Posterior Cervical Spine (PCS) Support	15
MAYFIELD® General Purpose Headrest	15

MAYFIELD® Infinity Support System	16
MAYFIELD® Infinity Support system	16
MAYFIELD® Infinity Coupler	16
MAYFIELD® Infinity Horsehoe	17
MAYFIELD® Infinity General Purpose Headrest	17
MAYFIELD® Infinity Skull Clamp System	18
MAYFIELD® Infinity Skull Clamp	18
MAYFIELD® Torque screw	18
MAYFIELD® Rocker Arms	19
MAYFIELD® Single Pin Accessory	19
BUDDE® Halo Retractor System	20
BUDDE® Halo Retractor System	20
Integra Brain retractor System	21
Integra Brain retractor System	21
MAYFIELD® Accessories	22
MAYFIELD® NeuroGen Adaptor	22
MAYFIELD® Crossbar Adaptor	22
MAYFIELD® Universal Side Rail Fitting	22
Integra™ Cranial Stabilisation Covers	22
MAYFIELD® Instrument Wall Storage Unit	23
MAYFIELD® Horseshoe Conversion Adaptor	23
MAYFIELD® Triad Horseshoe Conversion Adaptor with Knobs	23
BUDDE® Halo Side Rail Support	23
MAYFIELD® Infinity Radiolucent System	24
MAYFIELD® Infinity Radiolucent system	24
MAYFIELD® Removable Rocker Arm	24
MAYFIELD® Force Module mounted on a Metal-free Conversion accessory	24
MAYFIELD® Infinity XR2 Base Unit	25
MAYFIELD® Infinity XR2 Skull Clamp with improved swivel lock knob	25
MAYFIELD® Infinity XR2 Tri-Star Swivel Adaptor	25
MAYFIELD® Radiolucent System	26
MAYFIELD® Radiolucent Horseshoe Headrest	26
BUDDE® Halo Radiolucent Retractor System	27
BUDDE® Halo Radiolucent Retractor system	27
Index	28

Pricing:

Prices stated in specific written quotations are firm for thirty days from the date given, and are otherwise subject to change without prior notice. Pricing terms stated on written agreements are governed by such agreements.

Minimum Order Requirements:

Minimum order requirement is €250 (a €20 handling charge will be added to any order under €250).

Ordering Procedure:

A written purchase order on the customer's form may be requested for purchases of Integra LifeSciences products.

Acceptance of Orders:

Orders are accepted upon approval by Integra Customer Service.

Return Policy:

- Authorization from customer service must be obtained prior to returning product.
- Sterile product must be returned in unopened, undamaged cartons, packed to prevent damage.
- Non-sterile product must be returned in unused saleable condition in original package.
- Custom or special order products will not be accepted for credit.
- Credit will be issued for goods returned prior to ninety days from ship date with a 20% restocking charge. This assumes that the product returned is not damaged and can be verified to have not been used or opened.

Customer Service Contact Information:

International

+33 (0) 437 47 59 50
+33 (0) 437 47 59 25 (Fax)
csemea@integralife.com

France

+33 (0) 437 47 59 10
+33 (0) 437 47 59 29 (Fax)
custsvcfrence@integralife.com

Benelux

+32 (0) 2 257 4130
+32 (0) 2 253 2466 (Fax)
custsvcbenelux@integralife.com

United Kingdom

+44 (0) 1264 345 780
+44 (0) 1264 363 782 (Fax)
custsvcsuk@integralife.com

Germany

+49 (0) 2102 5535 6200
+49 (0) 2102 5536 636 (Fax)
custsvcgermany@integralife.com

Switzerland

+41 22 721 23 30
+41 22 721 23 99 (Fax)
custsvcsuisse@integralife.com

Integra Limited Warranty

INTEGRA LIFESCIENCES CORPORATION and its wholly owned subsidiaries ("INTEGRA") warrant to INTEGRA authorized distributors and the original purchaser only that each new INTEGRA product is free from manufacturing defects in material and workmanship under normal use and service for a period of one (1) year (except as otherwise expressly provided as to accessory items) from the date of delivery by INTEGRA (or its authorized distributor) to the original purchaser, but in no event beyond the expiration date stated on any product labeling. For purposes of products sold by INTEGRA through an authorized distributor of INTEGRA, "original purchaser" shall include the purchaser of INTEGRA products to whom the distributor first sells the product.

- Surgical instruments are guaranteed to be free from defects in material and workmanship when maintained and cleaned properly and used normally for their intended purpose.
- Any covered product that is placed by INTEGRA under a lease, rental or installment purchase agreement and that requires repair service during the term of such placement agreement shall be repaired in accordance with the terms of such agreement.

If any covered defect occurs during the warranty period or term of such placement agreement, the purchaser or distributor should communicate directly with INTEGRA. If purchaser or distributor seeks to invoke the terms of this warranty, the product must be returned to INTEGRA. The defective product should be returned promptly, properly packaged and postage prepaid. Loss or damage in return shipment to INTEGRA shall be at sender's risk. INTEGRA's sole responsibility under this warranty shall be repair or replacement, at INTEGRA's sole discretion at INTEGRA's expense, subject to the terms of this warranty and applicable agreements.

IN NO EVENT SHALL INTEGRA BE LIABLE FOR ANY INCIDENTAL, INDIRECT, CONSEQUENTIAL OR PUNITIVE DAMAGES IN CONNECTION WITH THE ACQUISITION OR USE OF ANY INTEGRA PRODUCT. Further, this warranty shall not apply to, and INTEGRA shall not be responsible for, any loss arising in connection with the purchase or use of any INTEGRA product that has been repaired by anyone other than an authorized INTEGRA service representative or altered in any way so as, in INTEGRA's judgment, to affect its stability or reliability, or which has been subject to misuse, negligence or accident, or which has been used otherwise than in accordance with the instructions furnished by INTEGRA.

THIS INTEGRA LIMITED WARRANTY IS EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, AND OF ALL OTHER OBLIGATIONS OR LIABILITIES ON INTEGRA'S PART OR THE PART OF ITS DISTRIBUTORS, AND INTEGRA NEITHER ASSUMES NOR AUTHORIZES ANY REPRESENTATIVE OR OTHER PERSON TO ASSUME FOR IT ANY OTHER LIABILITY IN CONNECTION WITH INTEGRA'S PRODUCTS.

INTEGRA DISCLAIMS ALL OTHER WARRANTIES, EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR OF FITNESS FOR A PARTICULAR PURPOSE OR APPLICATION OR WARRANTY OF QUALITY AS WELL AS ANY EXPRESS OR IMPLIED WARRANTY TO PATIENTS. No warranty or guarantee may be created by any act or statement nor may this Standard Warranty be modified in any way, except as a result of a writing signed by an officer of INTEGRA. These limitations on the creation or modification of this warranty may not be waived or modified orally or by any conduct.

IN NO EVENT SHALL INTEGRA AUTHORIZED DISTRIBUTORS BE LIABLE TOWARDS THE ORIGINAL PURCHASER FOR ANY INCIDENTAL, INDIRECT, CONSEQUENTIAL OR PUNITIVE DAMAGES IN CONNECTION WITH THE ACQUISITION OR USE OF ANY INTEGRA PRODUCT. Further, this warranty shall not apply to, and INTEGRA authorized distributors shall not be responsible towards the original purchaser for, any loss, arising in connection with the purchase or use of any INTEGRA product that has been repaired by anyone other than an authorized INTEGRA service representative or altered in any way so as to affect its stability or reliability, or which has been subject to misuse, negligence or accident, or which has been used otherwise than in accordance with the instructions furnished by INTEGRA.

THIS INTEGRA DISTRIBUTOR LIMITED WARRANTY IS EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES TOWARDS THE ORIGINAL PURCHASER, EXPRESS OR IMPLIED, AND OF ALL OTHER OBLIGATIONS OR LIABILITIES TOWARDS THE ORIGINAL PURCHASER ON INTEGRA AUTHORIZED DISTRIBUTOR'S PART. INTEGRA AUTHORIZED DISTRIBUTORS DISCLAIM ALL OTHER WARRANTIES TOWARDS THE ORIGINAL PURCHASER, EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR OF FITNESS FOR A PARTICULAR PURPOSE OR APPLICATION OR WARRANTY OF QUALITY AS WELL AS ANY EXPRESS OR IMPLIED WARRANTY TO PATIENTS.

MAYFIELD® surgical devices A heritage of innovation and excellence



Developed by neurosurgeon Dr. Frank Mayfield; in 1968, our MAYFIELD® Headrest System is the gold standard in cranial stabilization equipment.



MAYFIELD® Triad Skull Clamp
A1108

MAYFIELD®, the one headrest worthy of the name.

- Over 40 years of excellent service.
- Continual product development and improvement.
- New, innovative products.

Illustration of a complete MAYFIELD® Skull Clamp System

MAYFIELD® Modified Skull Clamp (A1059) shown with the Ultra Base Unit (A2101), Swivel Adaptor (A1018), and Adult Disposable Skull Pins (A1047).



MAYFIELD® Ultra 360 Patient Positioning System

New Support/Mounting device that facilitates attachment of MAYFIELD® Skull Clamps and Headrest System (A1012) to the OR table.

- Independent rotating and self-locking handles provide more fluid and flexible movement allowing better connection to the desired skull clamp position.
- Easier, more fluid movement for quicker positioning with skull clamp.
- Two double-cam action locking handles for easier opening and closing for secure, quick fixation.
- Set-up time is reduced because of the quick action handles.
- More flexible movement with added rotation for additional degrees of freedom.
- Unit can be placed in the reverse position to allow more reach and added height.
- Both a standard 360 Swivel Adaptor and a 360 Tri-star Swivel Adaptor are included.

Catalog No. A2009

System includes:

Component	Catalog #	Quantity
360 Adaptor Assembly	41A1614	1
360 Tri-Star Adaptor Assembly	41B1615	1



41A1614



41B1615

MAYFIELD® Ultra Base Unit

Support/Mounting device that facilitates attachment of MAYFIELD® Skull Clamps and Headrest Systems to the OR table.

- 33% more weight and force bearing capacity for improved patient safety.

Catalog No. A2101



MAYFIELD® Swivel Adaptor

Attaches to MAYFIELD® Base Unit to provide optimal flexibility in positioning via 360-degree rotation.

- Adapts headrests and skull clamps to MAYFIELD® Base Units.

Catalog No. A1018




MAYFIELD® Modified Swivel Adaptor

Same as A1018, with added spring-loaded pin; contacts ratchet arm of skull clamp to minimize all movement of the skull clamp components.

- Ideal for use with Image Guided Surgery (IGS) procedures.
- Compatible with all MAYFIELD® Skull Clamp, except the A1114 Infinity Skull Clamp.

Catalog No. A1115

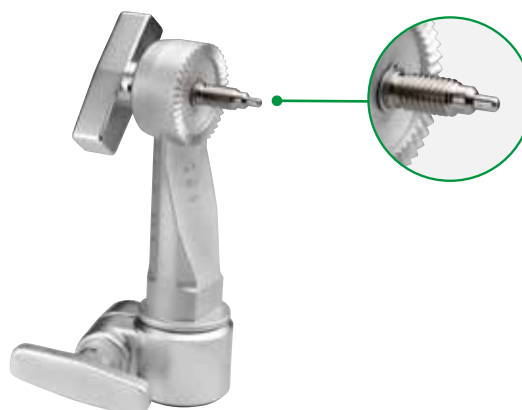

MAYFIELD® Ball Socket Swivel Adaptor

Attaches to MAYFIELD® Base Unit to provide optimal flexibility in positioning via 360-degree rotation with tilt in any axis.

- Unique locking pivot permits precise patient positioning.
- Exclusive spring-loaded mechanism provides assured stability and secure lock.
- Telescoping handle for easy access for tightening and loosening.

Catalog No. A1064





MAYFIELD® Tri-Star Swivel Adaptor

Attaches to MAYFIELD® Ultra Base Unit (A-2101) to provide optimal flexibility in positioning via 360-degree rotation. Tri-Star must be used with the Ultra Base Unit, the locking levers on both are gold-coloured for easy and convenient matching.

- Provides two additional starbursts to facilitate attachment of ancillary components such as Image Guided Surgery (IGS) dynamic reference frames and articulating arms.
- Starburst can rotate 360-degrees independent of skull clamp.


Catalog No. A2008

MAYFIELD® Modified Tri-Star Swivel Adaptor

Same as A2008, with added spring-loaded pin; contacts ratchet arm of skull clamp to minimize all movement of the skull clamp components.

- Compatible with all MAYFIELD® skull clamps, except the A1114 Infinity Skull Clamp.

Catalog No. A1116

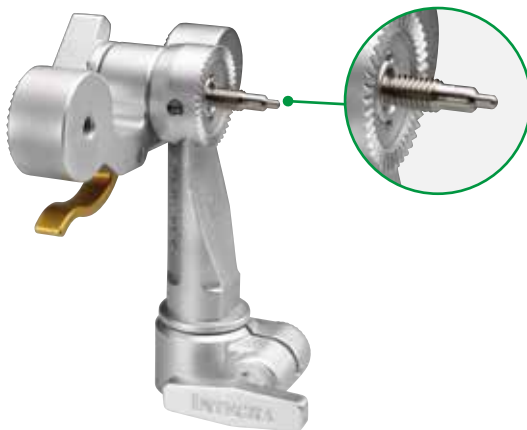

MAYFIELD® Infinity Skull Clamp System

Provides 3-point rigid cranial fixation.

- Added length and curved contours of the skull clamp uprights allow clamp to better interface with the horseshoe headrests of the Infinity Support System.
- Multiple rocker arm options:
Standard, Triad, Paediatric, Single-pin.
- 80 lb and 18 lb Torque Screw options.

Catalog No. A1114



See p.18 for complete system description.

MAYFIELD® Triad Skull Clamp

Provides 3-point rigid cranial fixation.

- Provides visual indication of clamping force on ALL THREE skull pins.
- Ability to adjust the force on ALL THREE skull pins.
- Equal impingement of pins ensures firm skull positioning and fixation.

Catalog No. A1108



MAYFIELD® Modified Skull Clamp

Provides 3-point rigid cranial fixation.

- Flexibility and safety in positioning of skull pins around critical areas.
- Quick-release lock.

Catalog No. A1059

**MAYFIELD® 2000 Skull Clamp**

Provides 3-point rigid cranial fixation.

- Streamlined body design and low profile single pin side facilitates greater exposure to the operative site.
- Unique swivel rocker arm mechanism permits greater flexibility in patient positioning.

Catalog No. A2000



MAYFIELD® Adult Reusable Skull Pins (Steel)

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Black O-ring connotes a reusable product; ensures a secure fit.

Catalog No. A1047



MAYFIELD® Child Reusable Skull Pins (Steel)

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Pinpoint reduced to control penetration when used with a child.
- Black O-ring connotes reusable product; ensures a secure fit.

Catalog No. A1048



MAYFIELD® Adult Disposable Skull Pins (Steel)

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Red O-ring connotes a disposable product; ensures a secure fit.

Catalog No. A1083



MAYFIELD® Child Disposable Skull Pins (Steel)

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Pinpoint reduced to control penetration when used with a child.
- Red O-ring connotes a disposable product; ensures a secure fit.

Catalog No. A1084



MAYFIELD® Adult Disposable Skull Pins

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- “Wings” design of the pin body improves the fit into receptacle holes of all MAYFIELD® Skull Clamps.
- Unique finger grooves molded into the pin body facilitate sure grip and easy removal from pin holders.
- Plastic base steel tip.

Catalog No. A1072





MAYFIELD® Adult Disposable Radiolucent Skull Pins

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Engineered from synthetic Sapphire material to significantly limit artifact during CT and MR imaging.
- “Wings” design of the pin body improves the fit into receptacle holes of all MAYFIELD® Skull Clamps.
- Unique finger grooves molded into the pin body facilitate sure grip and easy removal from pin holders.

Catalog No. A2020

NEW



MAYFIELD® Adult Disposable Titanium Skull Pins

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Designed and tested for use in MR systems up to 3.0 Tesla (T)*.
- “Wings” design of the pin body improves the fit into receptacle holes of all MAYFIELD® Skull Clamps.
- Unique finger grooves molded into the pin body facilitate sure grip and easy removal from pin holders.

Catalog No. A1120

NEW



MAYFIELD® Child Disposable Titanium Skull Pins

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Pinpoint reduced to control penetration when used with a child.
- Designed and tested for use in MR systems up to 3.0 Tesla (T)*.
- “Wings” design of the pin body improves the fit into receptacle holes of all MAYFIELD® Skull Clamps.
- Unique finger grooves molded into the pin body facilitate sure grip and easy removal from pin holders.

Catalog No. A1119

NEW



MAYFIELD® Adult Reusable Titanium Skull Pins

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Black O-ring connotes a reusable product; ensures a secure fit.
- Designed and tested for use in MR systems up to 3.0 Tesla (T)*.

Catalog No. A1121

NEW



MAYFIELD® Child Reusable Titanium Skull Pins

Fit all MAYFIELD® Skull Clamps to allow rigid cranial fixation.

- Pinpoint reduced to control penetration when used with a child.
- Designed and tested for use in MR systems up to 3.0 Tesla (T)*.
- Black O-ring connotes a reusable product; ensures a secure fit.

Catalog No. A1122

Illustration of a complete MAYFIELD® Headrest System



MAYFIELD® Adult Horseshoe Headrest (A1012) with Ball Socket Swivel Adaptor (A1064) and Ultra Base Unit (A2101).

MAYFIELD® Adult Horseshoe Headrests

Horseshoe shaped headrest provides cranial support in the prone or supine position for adults and older children.

- Vertical and lateral adjustments allow flexibility in patient positioning.
- Adjustable pad base accommodates various head sizes.
- Fluid resistant gel pads provide comfort and reduce the incidence of pressure necrosis.
- Pulley Rod attachment for skeletal traction.
- Can be attached to the transitionnal member of the MAYFIELD® Ultra Base unit A2101 thanks to the A1110 MAYFIELD® Triad Horseshoe conversion adaptor with knobs (see page 19).

Catalog No. A1012



A1012

- Same as A1012 except connection is made directly to the Transitional member for chin-to-chest positioning.
- Can be attached to the MAYFIELD® swivel adaptors or A2009 MAYFIELD® Ultra 360 patient positioning system thanks to the A1109 MAYFIELD® Horseshoe conversion adaptor (see page 19 of the catalogue).

Catalog No. A1011



A1011

MAYFIELD® Paediatric Horseshoe Headrest

Horseshoe shaped headrest provides cranial support in the prone or supine position for infants and young children.

- Vertical and lateral adjustments allow flexibility in patient positioning.
- Adjustable pad base accommodates various head sizes.
- Fluid resistant gel pads provide comfort and reduce the incidence of pressure necrosis.
- Easily cleaned.
- Can be attached to the MAYFIELD® swivel adaptors or A2009 MAYFIELD® Ultra 360™ patient positioning system thanks to the A1109 MAYFIELD® Horseshoe conversion adaptor (see page 19 of the catalogue).

Catalog No. A1051



Papadopoulos Posterior Cervical Spine (PCS) Support

Supports the cervical spine in an anatomical position.

- Firm posterior cervical spine support for anterior approaches.
- Ratcheting System allows for easy intraoperative cervical extension.
- Base is radiolucent to facilitate intraoperative imaging.

Catalog No. A1073



MAYFIELD® General Purpose Headrest

Provides cranial support and positioning when rigid fixation is not needed.

- Headrest pads easily contour to a patient's head.

Catalog No. A1008



MAYFIELD® Infinity Support System

Provides multiple headrests for cranial support in the prone or supine position for paediatric patients (ages 5+) and adults.

The MAYFIELD® Infinity attaches to any MAYFIELD® Base Unit, and creates a method for simultaneous application of a head support (adult, paediatric or general purpose) and a skull clamp. This is especially important for patients who require a low clamping force to stabilize the position of the head, such as paediatric patients (ages 5+) and trauma patients.

The Infinity head supports can be positioned close to the O.R. Table, which is important for patients with short neck lengths.

Catalog No. A1112



Infinity system with the Ultra Base Unit A2101.

System includes:

Component	Catalog #	Quantity
Infinity Support Coupler	419B1010	1
Horseshoe – Paediatric	419A1030	1
Horseshoe – Adult	419A1031	1
Horseshoe – General Purpose	419A1050	1
Mini Gel Pad for General Purpose	419A1048	2
Storage Case (not pictured)	419A1060	1

Replacement Gel Pads:

Component	Catalog #	Quantity
Paediatric – Left	440C1011	1
Paediatric – Right	440C1012	1
Adult – Left	41C1444	1
Adult – Right	41C1445	1
Mini Gel Pad for General Purpose	419A1048	1

Illustration of the Infinity Support System mounted on the Ultra Base Unit A2101, the MAYFIELD® Swivel Adaptor A1018 and the MAYFIELD® Triad Skull Clamp A1108.



MAYFIELD® Infinity Coupler (not sold separately)

- Coupler easily mounts to any MAYFIELD® Base Unit.
- Head supports can be adjusted vertically and laterally for flexibility in patient positioning, then securely locked in place.
- Position of head support can be tilted on ball socket joint and locked in place.

Catalog No. 419B1010



**MAYFIELD® Infinity Adult Horseshoe (not sold separately)**

- Provides cranial support for adults and older children in either prone or supine positions.
- Adjustable pad base (faceplates) accommodates various head sizes.
- Fluid resistant gel pads provided for comfort and reduction in the incidence of pressure necrosis.
- Easily cleaned.

Catalog No. 419A1031

**MAYFIELD® Infinity Paediatric Horseshoe (not sold separately)**

- Provides cranial support in either prone or supine positions for infants and young children.
- Same benefits of the MAYFIELD® Infinity Adult Horseshoe (419A1031).

Catalog No. 419A1030

**MAYFIELD® Infinity General Purpose Headrest (not sold separately)**

- Provides cranial support in either supine or semi-sitting procedures.
- Mini gel pads can be positioned in multiple receptacles and can be tilted on their hinges.

Catalog No. 419A1050



MAYFIELD® Infinity Skull Clamp

Multiple, quick-change, snap-in rocker arm options provide the ability to match the stabilization points to the size of the patient's head for optimal fixation and surgical access. No tools are needed to change rocker arms.

Catalog No. A1114



MAYFIELD® Infinity Skull Clamp System (A1114) includes:

Component	Catalog #	Quantity
Infinity Skull Clamp	A1114	1 each
Adult Rocker Arm	419B1204	1 each
Child's Rocker Arm	419B1205	1 each
MAYFIELD® Triad Rocker Arm	419B1206	1 each
Single Pin Accessory	419A1207	1 each
Reduced Load 18-lb o	419A1208	1 each
80-lb Torque Screw	419A1209	1 each
Adjustment Wrench	41A1482	1 each
Storage Case	419A1250	1 each

Infinity Skull Clamp shown with:

- 80-lb Torque Screw
- Standard Adult Rocker Arm
- Infinity Support System A1112
- MAYFIELD® Ultra Base Unit A2101
- MAYFIELD® Tri-Star Swivel Adaptor A2008

MAYFIELD® Torque Screw

Along with the standard MAYFIELD® 80-lb Torque Screw, a reduced load 18-lb Torque Screw is provided for procedures where lower, more controlled pressure and less penetration of the pins is desired.

(Note: The reduced load 18-lb Torque Screw should always be used in conjunction with a horseshoe headrest or other form of primary patient head support.)



80-lb Torque Screw
Catalog No. 419A1209



18-lb Torque Screw
Catalog No. 419A1208



MAYFIELD® Triad Rocker Arm

Provides visual indication of clamping force of the two pins on the rocker arm. Also has the ability to adjust clamping force on each pin. Set-up is repeatable, making training of the application of skull clamps relatively easy.

Catalog No. 419B1206



MAYFIELD® Standard Adult Rocker Arm

The same two-pin rocker that is provided with the standard MAYFIELD® skull clamp.

Catalog No. 419B1204



MAYFIELD® Child's Rocker Arm

Rocker Arm configured for smaller size patient skulls.

Catalog No. 419B1205



MAYFIELD® Single Pin Accessory

Provides a directly opposing point of stabilization with the pin in the Torque Screw. Flexibility in positioning of pins around critical areas and the surgical site.

Catalog No. 419A1207

BUDDE® Halo Retractor System

Attaches to MAYFIELD® Skull Clamps or to the side rails of the OR table to provide a sturdy platform for retraction or hand rest for the surgeon.

- Double-Hinged ContinuTrac halo ring provides 360-degrees of uninterrupted retractor arm placement for ultimate flexibility.
- Low profile design allows for ease of instrument passage within the surgical field.
- Ideal for posterior fossa procedures.
- Blades can be micro-manipulated without flex arm extension readjustment.

Catalog No. A1040



Shown on a MAYFIELD® Triad Skull Clamp (A1108, see on page 5) and A2008 MAYFIELD® Tri-Star Adaptor Modified (seen on page 9).

System includes:

Component	Catalog #	Quantity
ContinuTrac™ Halo Ring	438B1178	1
9" (229 mm) Halo Flex Arm	438A1011	3
BUDDE® Halo Support Bracket	438B1010	2
BUDDE® Halo Support Rod (5/16")	438A1188	2
BUDDE® Halo Pattie Tray	438A1037	1
BUDDE® Halo Mini Vise - SS (Not Shown)	438A1140	1
BUDDE® Curved Retractor	438A1077	2
TEW Micro Retractor Blades - SS	438A1082	1
TEW Micro Retractor Blades - SS	438A1083	1
Halo Adjustment Wrench	438A1024	2
Halo Poly Vac Sterilization Case	438A1096	1

Optional Accessories:

Component	Catalog #
12" (300 mm) Halo Flex Arm - SS	438A1012
6" (150 mm) Halo Flex Arm - SS	438A1160
Hooked Micro Retractor Blades (2 mm) - SS	438A1151
Hooked Micro Retractor Blades (4 mm) - SS	438A1152
BUDDE® Tapered Retractor Blades - SS	438A1154
Serrated Micro Retractor Blades (2 mm) - SS	438A1155
Serrated Micro Retractor Blades (4 mm) - SS	438A1157
Samy Vascular Retractors Short (3, 5, 7 mm)	438A1200A
Samy Vascular Retractors Long (3, 5, 7 mm)	438A1200B



Integra Brain retractor system

The next generation in Retractor possibilities. This is a complete system with individual components providing as much or as little retraction as needed for each procedure.

- Floating Secondary Bars provide the versatility of the ball-socket for positioning the bar where it is needed.
- Only a Primary Bar, Floating Secondary Bar, Retractor Arm, and a Retractor Blade is needed for retraction, providing a precious ease of set up.

Catalog No. A2012



System includes: SIDE RAIL ATTACHMENT

- For those cases where a skull clamp is not utilized, the Integra Brain Retractor System can be attached directly to the side rails of the OR table
- By connecting the Integra Brain Retractor System to the side rails of the OR table, the Integra Brain Retractor System can also be used for spine cases as well

Component	Catalog #	Quantity
Primary bars and clamp assembly	417A1000	2 each
Floating Secondary Bars	417A1001	4 each
Retractor Arms 9.5"	417A1002	4 each
Handrest	417A1003	1 each
PATTIE Tray	417A1004	1 each
Tapered Retractor Blade, 1/4" to 3/4"(6.4 mm to 19.1 mm)	417A1005	1 each
Tapered Retractor Blade, 5/32" to 5/8"(4.0 mm to 15.9 mm)	417A1006	1 each
Tapered Retractor Blade, 1/8" to 1/2"(3.2 mm to 12.7 mm)	417A1007	1 each
Tapered Retractor Blade, 3/32" to 3/8"(2.4 mm to 9.5 mm)	417A1008	1 each
Tapered Retractor Blade, 1/16" to 1/4"(1.6 mm to 6.4 mm)	417A1009	1 each
Sterilization Case	417A1010	1 each

MAYFIELD® NeuroGen Adaptor

Support/Mounting device that allows attachment of base units to the OR table when receptacle holes are absent or inappropriate.

- Fits a wide variety of table widths.
- Easily attaches to OR table side rails.
- Quick, adjustable end bracket to fit a wide range of table widths.
- Provides greater lateral adjustment with multiple mounting holes for base unit attachment.
- Several holes to better laterally adjust the set up of the system.

Catalog No. A1113



MAYFIELD® Crossbar Adaptor

- For seated position procedures, easily mounts to the OR table side rails to adapt all MAYFIELD® Skull Clamps for posterior fossa craniotomy procedures or posterior cervical laminectomies.
- Allows for unobstructed view of the posterior cervical spine and occipital areas.
- Compatible with most currently available OR tables.


Catalog No. A1015



MAYFIELD® Universal Side Rail Fitting

- For attaching the MAYFIELD® Crossbar Adaptor (A1015) and other products to the OR table.
- Compatible with most currently available OR tables.

Catalog No. A1060

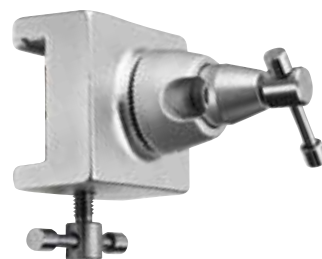

INTEGRA™ Cranial Stabilisation Covers

Non-Sterile, Disposable Covers designed to protect MAYFIELD® Systems.

- Reduce cleaning time.
- Protect critical components.
- Provide patients and staff an increased level of safety.

Catalog No. A1118





MAYFIELD® Instrument Wall Storage Unit

Wall board used to organize and store MAYFIELD® products.

- Printed silhouettes allow quick and easy storage.
- Holds all MAYFIELD® Headrest components.

Catalog No. A1023



MAYFIELD® Horseshoe Conversion Adaptor

Allows an A1011 MAYFIELD® Horseshoe Headrest and A1051 MAYFIELD® Paediatric Horseshoe Headrest to attach to A1018, A1064, A2008 and Adaptor assembly of the A2009 MAYFIELD® Ultra 360™ Patient positioning system.

Catalog No. A1109



MAYFIELD® Triad Horseshoe Conversion Adaptor with Knobs

Allows an A1012 MAYFIELD® Horseshoe Headrest to attach to standard Transition Members.

- Provides a means of obtaining chin-to-chest positioning.

Catalog No. A1110



BUDDE® Halo Side Rail Support

Allows the BUDDE® Halo Retractor System to attach directly to the operating table.

Catalog No. A1067

MAYFIELD® Infinity Radiolucent System

MAYFIELD®, in a constant technological evolution strategy and always striving to bring new solutions to surgeons, re-design the MAYFIELD® Radiolucent product line.

Preservation of the target is key for visualization during imaging. Using metal cranial stabilization equipment during such surgeries results in unwanted artifacts in the image which hinder the capabilities and usefulness of images.

The new MAYFIELD® Infinity XR2 is part of the growing “Infinity” brand. Designed to upgrade the existing top-of-the-line Radiolucent Systems. All of the new components are backwards compatible to the existing Skull Clamp and Base Unit.



Total XR2 System:
- MAYFIELD® XR2 Skull Clamp A2114
- MAYFIELD® XR2 Swivel Adaptor
- MAYFIELD® XR2 Base Unit A2079

Removable Rocker Arm

- No tool is needed to remove the Rocker Arm for cleaning or changing to another size.
- Tapered shaft hold the Rocker Arm in place.

Catalog No. 439A1010



Removable Child Rocker Arm

- Rocker Arm holds two skull pins, allows for positioning on smaller skull sizes.

Catalog No. 439A1091



Force Module mounted on a Metal-Free Conversion Accessory

- Converts the skull clamp to be free of any metal in the product during imaging.

Catalog No. 439A1092 (metal-free conversion accessory)

Catalog No. 439A1093 (force module)





MAYFIELD® Infinity XR2 Base Unit

Improved Rigidity for a more stable operating platform for the surgeon.

- Set-up and positioning is easy to set-up with the simplified linkage & components saving OR and staff time.
- Enhanced level of rigid fixation with a wide imaging window for procedures using to better accommodate intra-operative imaging modalities.
- Colour Code : the Blue and Yellow Base Unit Screws make it easy to install and take down for cleaning.
- Includes a standard Radiolucent Swivel Adaptor

Catalog No. A2079

MAYFIELD® Infinity XR2 Skull Clamp with improved swivel lock knob



- 80-lbs Torque Knob on the single-pin side for a familiar position & easy to read.
- 1/4 turn Rocker Lock in a familiar style and easy to use.
- Removable Adult Rocker Arm for easy cleaning.
- Double-Pawl ratchet mechanism for easier application of the skull clamp by the surgeon to the patient.
- Includes the adult Rocker arm 439A1010

Catalog No. A2114



MAYFIELD® Infinity XR2 Tri-Star Swivel Adaptor

- Designed to provide flexibility in positioning the patient's head by allowing rotation in 360 degrees.*
- Provides two extra starburst attachments for ancillary equipment, when needed.**
- Allows for a strong attachment of image-guided surgery components.**
- Easy set-up and firm fixation onto the Base Unit and the Skull Clamp.*

Catalog No. A2111

* for both adaptors **only for A2111

MAYFIELD® Adult Radiolucent Horseshoe Headrest

- Can be used in conjunction with the MAYFIELD® Radiolucent System.
- Provides support in the prone and supine position during diagnostic examination or surgical procedures.
- Offers the same benefits of the standard Horseshoe Headrest (A1012).
- Composite carbon fiber provides strength, durability and minimal artifact during imaging.
- Designed to be used with the MAYFIELD® Infinity XR2 Base Unit A2079.

Catalog No. A2010



MAYFIELD® Paediatric Radiolucent Horseshoe Headrest

- Can be used in conjunction with the MAYFIELD® Radiolucent System.
- Provides support in the prone and supine position for children during diagnostic examination or surgical procedures.
- Offers the same benefits of the standard Horseshoe Headrest (A1012).
- Composite carbon fiber provides strength, durability and minimal artifact during imaging.

Catalog No. A2011





BUDDE® Halo Radiolucent A1096
shown on a XR2 Infinity Skull Clamp A2114.

BUDDE® Halo Radiolucent Retractor System

Attaches to MAYFIELD® Radiolucent Skull Clamps. Provides a sturdy platform for retraction and a ready support or hand rest for the surgeon.

- Radiolucent components allow the system to remain in place during intraoperative imaging for optimal efficiency.
- Compatible with the MAYFIELD® Radiolucent Headrest System (A2114), the BUDDE® Halo Radiolucent Retractor System is designed to significantly reduce or limit artefact during intracranial procedures.
- Imaging compatible titanium retractor arms can remain in place during scanning.
- Ring can be quickly changed to a half ring set-up for improved access.
- Blades can be micro-manipulated without flex arm extension readjustment.

Catalog No. A1096

System includes:

Component	Catalog #	Quantity
Halo ring	443A1000	1 chacun
9" (229 mm) halo flex arm - TI	443A1070A	3 chacun
Rod support clamp assembly	443A1060	2 chacun
Radiolucent support rod	443B1052	2 chacun
BUDDE® Halo Pattie tray - SS	438A1037	1 chacun
BUDDE® Halo Mini Vise - SS	438A1140	1 chacun
BUDDE® Halo adjustment wrench	443A1004	1 chacun
BUDDE® Halo Poly Vac sterilization case	443A1005	1 chacun
Retractor arm support	443A1082	3 chacun
Samy vascular retractors short (3, 5, 7 mm)	438A1200A	3 chacun
Samy vascular retractors long (3, 5, 7 mm)	438A1200B	3 chacun

Optional Accessories:

Component	Catalog #
12" (300 mm) Halo flex arm - TI	443A1070B
6" (150 mm) Halo flex arm - TI	443A1070C
TEW hooked micro retractor blades (2 mm) - SS	438A1151
TEW hooked micro retractor blades (4 mm) - SS	438A1152
BUDDE® Tapered retractor blades - SS	438A1154
TEW serrated micro retractor blades (2 mm) - SS	438A1155
TEW serrated micro retractor blades (4 mm) - SS	438A1157

Products
References

#

360 Adaptor Assembly.....	7
360 Tri-Star Adaptor Assembly.....	7
2000 Skull Clamp.....	11

A

Accessories.....	22, 23
Adjustment Wrench.....	18
Adult Disposable Skull Pins.....	12
Adult Disposable Titanium Skull Pins.....	13
Adult Horseshoe Headrest.....	14
Adult Radiolucent Horseshoe Headrest.....	26
Adult Reusable Skull Pins.....	12
Adult Reusable Titanium Skull Pins.....	13
Adult Rocker Arm.....	18

B

Ball Socket Swivel Adaptor.....	8, 14
Base Units.....	6, 7
BUDDE® Curved Retractor.....	20
BUDDE® Halo.....	27
BUDDE® halo adjustment wrench.....	27
BUDDE® halo mini vise.....	27
BUDDE®Halo Mini Vise.....	20
BUDDE® halo pattie tray.....	27
BUDDE®Halo Pattie Tray.....	20
BUDDE® Halo Radiolucent.....	27
BUDDE® Halo Radiolucent Retractor System.....	20, 27
BUDDE® Halo Retractor System.....	20
BUDDE® Halo Side Rail Support.....	23
BUDDE® Halo Support Bracket.....	20
BUDDE®Halo Support Rod.....	20
budde® tapered retractor blades.....	27
BUDDE® Tapered Retractor Blades.....	20

C

Child Disposable Skull Pins (Steel).....	12
Child Disposable Titanium Skull Pins.....	13
Child Reusable Skull Pins.....	12
Child Reusable Titanium Skull Pins.....	13
Child's Rocker Arm.....	18, 19
ContinuTrac™ Halo Ring.....	20
Cranial Stabilisation Covers.....	22
Crossbar Adaptor.....	22

D

Disposable Skull Pins.....	6
----------------------------	---

F

Floating Secondary Bars.....	21
Force Module.....	24

G

Gel Pad.....	16
General Headrests.....	15
General Purpose Headrest.....	15

H

Halo Adjustment Wrench.....	20
Halo Flex Arm.....	20, 27
Halo Poly Vac sterilization case.....	27
Halo Poly Vac Sterilization Case.....	20
Halo ring.....	27
Handrest.....	21
Headrest System.....	14
Horseshoe – Adult.....	16
Horseshoe Conversion Adaptor.....	23
Horseshoe – General Purpose.....	16
Horseshoe Headrests.....	14
Horseshoe – Paediatric.....	16

I

Infinity Adult Horseshoe.....	17
Infinity Coupler.....	16
Infinity General Purpose Headrest.....	17
Infinity Paediatric Horseshoe.....	17
Infinity Radiolucent System.....	24, 25
Infinity Skull Clamp.....	10, 18
Infinity Skull Clamp System.....	18, 19
Infinity Support Coupler.....	16
Infinity Support System.....	16, 17, 18
Infinity XR2 Base Unit.....	25
Infinity XR2 Skull Clamp with improved swivel lock knob.....	25
Infinity XR2 Tri-Star Swivel Adaptor.....	25
Instrument Wall Storage Unit.....	23
Integra Brain retractor system.....	21

M

MAYFIELD® Triad.....	5, 10, 16, 18, 19, 20, 23
Metal-Free Conversion Accessory.....	24
micro retractor blade.....	27
Micro Retractor Blade.....	20

N

NeuroGen Adaptor.....	22
-----------------------	----

P

Paediatric Horseshoe Headrest	15
Paediatric Radiolucent Horseshoe Headrest	26
Papadopoulos Posterior Cervical Spine (PCS) Support.....	15
PATTIE Tray	21
Primary Bars.....	21

R

Radiolucent Retractor System	27
Radiolucent Skull Pins	13
Radiolucent Support Rod	27
Radiolucent System.....	26
Removable Child Rocker Arm	24
Removable Rocker Arm.....	24
Retractor Arms	21
Retractor Arm Support.....	27
Rod Support Clamp Assembly	27

S

Samy vascular retractors	27
Single Pin Accessory	18, 19
Skull Clamp.....	6, 10, 11
Skull Clamp System.....	6
Skull Pins.....	12, 13
Standard Adult Rocker Arm.....	18, 19
Sterilization Case	21
Storage Case.....	16, 18
Swivel Adaptor	6, 8, 9, 16

T

Tapered Retractor Blade.....	21
TEW™ Micro Retractor Blades	20
Torque Screw	18
Tri-Star Adaptor	20
Tri-Star Swivel Adaptor	18

U

Ultra 360 Patient Positioning System.....	7
Ultra Base Unit	6, 7, 14, 16, 18
Universal Side Rail Fitting.....	22

V

Vascular Retractors.....	20
--------------------------	----

X

XR2 Base Unit.....	24
XR2 Skull Clamp	24
XR2 Swivel Adaptor.....	24
XR2 System.....	24

41A1614	7	443A1004	27
41A1615	7	443A1005	27
41B1615	7	443A1060	27
41C1444	16	443A1070A	27
41C1445	16	443A1070B	27
417A1000	21	443A1070C	27
417A1001	21	443A1082	27
417A1002	21	443B1052	27
417A1003	21		
417A1004	21	A1008	15
417A1005	21	A1011	14
417A1006	21	A1012	7, 14
417A1007	21	A1015	22
417A1008	21	A1018	6, 8, 16
417A1009	21	A1023	23
417A1010	21	A1040	20
419A1030	16, 17	A1047	6, 12
419A1031	16, 17	A1048	12
419A1048	16	A1051	15
419A1050	16, 17	A1059	6, 11
419A1060	16	A1060	22
419A1207	19	A1064	8, 14
419A1208	18	A1067	23
419A1209	18	A1072	12
419B1010	16	A1073	15
419B1204	19	A1083	12
419B1205	19	A1084	12
419B1206	19	A1096	27
438A1011	20	A1108	5, 10, 16
438A1012	20	A1109	23
438A1024	20	A1110	23
438A1037	20, 27	A1112	16, 18
438A1077	20	A1113	22
438A1082	20	A1114	10, 18
438A1083	20	A1115	8
438A1096	20	A1116	9
438A1140	20, 27	A1118	22
438A1151	20, 27	A1119	13
438A1152	20, 27	A1120	13
438A1154	20, 27	A1121	13
438A1155	20, 27	A1122	13
438A1157	20, 27	A2000	11
438A1160	20	A2008	9, 18
438A1188	20	A2009	7
438A1200A	20, 27	A2010	26
438A1200B	27	A2011	26
438B1010	20	A2012	21
438B1178	20	A2020	13
439A1010	24	A2079	24, 25
439A1091	24	A2101	6, 7, 14, 16, 18
439A1092	24	A2111	25
439A1093	24	A2114	24, 25
440C1011	16		
440C1012	16		
443A1000	27		



Distributed by

Integra LifeSciences Services (France) SAS
Sales & Marketing EMEA
Immeuble Séquoia 2 • 97 allée Alexandre Borodine
Parc technologique de la Porte des Alpes
69800 Saint Priest • FRANCE
+33 (0)4 37 47 59 00 phone • +33 (0)4 37 47 59 99 fax
emea.info@integralife.com • [integralife.eu](mailto:emea.info@integralife.com)

Customer Services

International: +33 (0)4 37 47 59 50 • +33 (0)4 37 47 59 25 (Fax) • csema@integralife.com
France: +33 (0)4 37 47 59 10 • +33 (0)4 37 47 59 29 (Fax) • custsvcfrence@integralife.com
United Kingdom: +44 (0)1 264 345 780 • +44 (0)1 264 363 782 (Fax) • custsvcs.uk@integralife.com
Germany: +49 (0) 2102 5535 6200 • +49 (0)2 102 5536 636 (Fax) • custsvcgermany@integralife.com
Benelux: +32 (0)2 257 4130 • +32 (0)2 253 2466 (Fax) • custsvcbenelux@integralife.com
Switzerland: +41 (0)2 27 21 23 30 • +41 (0)2 27 21 23 99 (Fax) • custsvcsuisse@integralife.com

Caution: Applicable laws restrict these products to sale by or on the order of a physician.

All the references numbers mentioned on this document are CE marked according to European council directive 93/42/EEC on medical devices, unless specifically identified as "NOT CE MARKED".

All text and diagrams contained in this catalog are intended as guidelines and are for illustrative purposes only. No information included in this catalog is intended to replace the Directions for Use packaged with each Integra product. Always refer to the appropriate Directions for Use for complete clinical instructions. For technical and additional clinical information on Integra products, please contact your Integra NeuroSpecialist.

Availability of these products might vary from a given country or region to another, as a result of specific local regulatory approval or clearance requirements for sale in such country or region.

BUDDE and the Integra logo are trademarks or registered trademarks of Integra LifeSciences Corporation or its subsidiaries. Microsensor and Codman are registered trademarks of their owners. MAYFIELD is a registered trademark of SM USA, Inc. and is used by Integra LifeSciences Corporation under license.